

# Caring for Your Catheter at Home

## A PATIENT'S GUIDE



**Chesapeake**  
**UROLOGY®**

*The most personal care for life's most personal issues.®*

**855-405-7100**  
**chesapeakeurology.com**

**Chesapeake**  
**UROLOGY®**

## This information will help you care for your urinary (Foley) catheter while you are at home.

You have had a Foley catheter (a thin, flexible tube) placed in your bladder to drain your urine. It is held inside your bladder by a balloon filled with water. The parts of the catheter outside your body are shown in Figure 1. Figure 2 illustrates the parts of the Foley catheter.


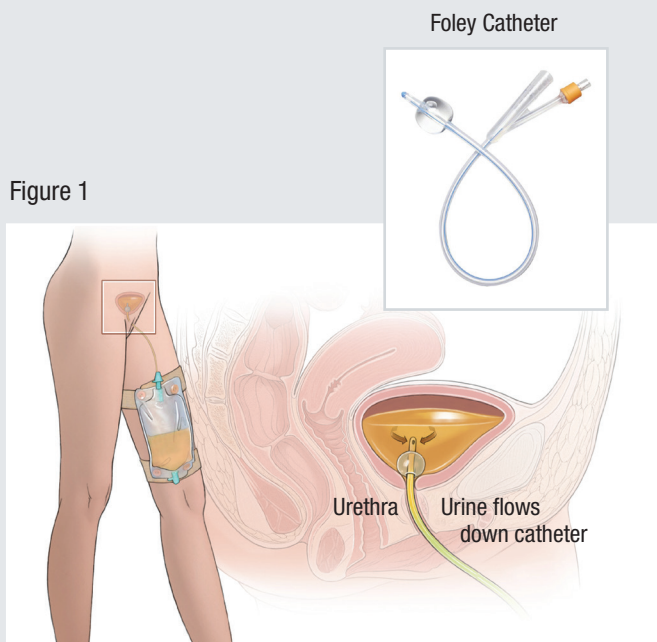


Figure 2

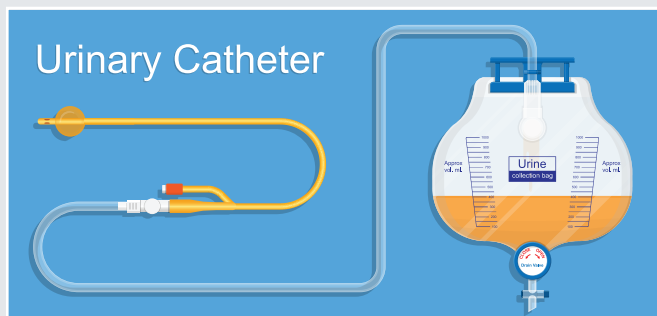

You need to perform catheter care every day. This involves cleaning the catheter, changing the drain collection bags, and washing the drainage bags. Follow these instructions to care for your catheter.

## CHANGING THE DRAINAGE BAG

You will need the following supplies:

- ☐ A clean washcloth (not one already used for bathing) or a 4"x 4" piece of gauze
- ☐ Night or leg drainage bag (whichever one you are switching to)
- ☐ 2 alcohol pads

1. Wash your hands thoroughly with soap and warm water.
2. Empty the urine from the drainage bag into the toilet. Make sure that the spout of the drainage bag never touches the side of the toilet or any emptying container.
3. Place the clean cloth or gauze under the connector to catch any leakage.
4. Pinch off the catheter with your fingers and disconnect the used bag.
5. Wipe the end of the catheter with an alcohol pad.
6. Wipe the connector on the new bag with the second alcohol pad.
7. Connect the clean bag to the catheter and release your finger pinch.
8. Check all connections. Straighten any kinks or twists in the tubing.

## CARING FOR THE DRAINAGE BAGS

### Caring for the leg bag

- Make sure to secure the leg bag to your leg with the elastic or velcro straps provided.
- If the straps leave a mark on your leg, they are too tight and should be loosened.
- Try to avoid letting the leg bag drop below the knee joint.
- Empty the leg bag into the toilet through the spout at the bottom every 2 to 4 hours, as needed.
- Don't let the bag become completely full.



### Caring for the night bag

- Always keep the night bag below the level of your bladder.
- To hang your night bag while you sleep, place a clean plastic bag inside of a wastebasket. Hang the night bag on the inside of the wastebasket.

### Cleaning the drainage bags

You will need the following supplies:

- ☐ White vinegar
- ☐ Cool water

1. Wash your hands thoroughly with soap and warm water.
2. Rinse the bag with cool water. Do not use hot water because it can damage the plastic equipment.
3. To decrease odor, fill the bag halfway with a mixture of 1 part white vinegar and 3 parts water. Shake the bag and let it sit for 15 minutes.
4. Rinse the bag with cool water and hang it up to dry.

## Preventing Infection

- Keep the drainage bag below the level of your bladder and off the floor at all times.
- Keep the catheter secured to your thigh to prevent it from moving.
- Don't lie on your catheter or block the flow of urine in the tubing.
- Shower daily to keep the catheter clean.
- Clean your hands before and after touching the catheter or bag.

## Things to Note

- You may see some blood or urine around where the catheter enters your body, especially when walking or having a bowel movement. This is normal, as long as there's urine draining into the drainage bag. If there's not, call your doctor or nurse.
- While you have your catheter, drink an extra glass of water with each meal and at bedtime. (Unless you take a "fluid pill" or diuretic medication.)

## Call Your Doctor Immediately if:

- Your catheter comes out. Do not try to replace it yourself.
- You have a temperature of 101 °F (38.3 °C) or higher.
- You are making less urine than usual.
- You have foul-smelling urine.
- You have bright red blood or large blood clots in your urine.
- You have abdominal (belly) pain and no urine in your catheter bag.

## Showering

You can shower while you have your Foley catheter. However, do not take a bath until after your catheter is removed. This is because taking a bath while you have the Foley catheter puts you at risk for infections.

Shower by removing the cloth leg straps from the drainage bag and letting the catheter bag gently hang while you wash.


## CLEANING YOUR CATHETER

You can clean your catheter while you are in the shower.

You will need the following supplies:

- ☐ Mild soap
- ☐ Water
- ☐ 1 Catheter leg strap

1. Wash your hands thoroughly with soap and water.
2. Using mild soap and water, clean your genital area.
  - Men should retract the foreskin, if needed, and clean the area, including the penis. Make sure that you replace the skin back in the normal position.
  - Women should separate the labia, and clean the area from front to back.
3. Clean your urethra (urinary opening), which is where the catheter enters your body.
4. Clean the catheter from where it enters your body and then down, away from your body. Hold the catheter at the point it enters your body so that you don't put tension on it.
5. Rinse the area well and dry it gently.
6. Use the leg strap provided to attach the catheter to your leg to keep it from moving.

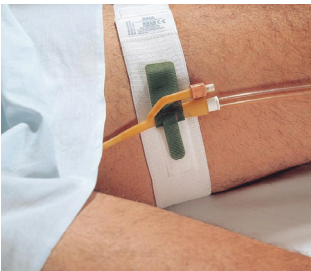

## HOW TO REMOVE YOUR CATHETER SAFELY:

1. Wash your hands thoroughly with soap and water.
2. Empty the drainage bag.
3. While sitting on the toilet, use scissors to cut the balloon port. This allows approximately a teaspoon of sterile water to automatically drain from the balloon inside of your bladder.
4. When the water is finished draining, gently slide the catheter out of your bladder. If you meet any resistance, please DO NOT pull with force. Please call your physician for direction.

